

ABHD8 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP00125
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from human ABHD8. AA range:291-340
Mol wt	47331
Species reactivity	Human, Mouse
Clonality	Polyclonal
Recommended application	WB, IHC-p, ELISA
Concentration	1 mg/ml
Full name	ABHD8 Antibody
Synonyms	ABHD8; Abhydrolase domain-containing protein 8

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

ABHD8 (abhydrolase domain containing 8) is upstream of, and in a head-to-head orientation with the gene for the mitochondrial ribosomal protein L34. The predicted protein contains alpha/beta hydrolase fold and secretory lipase domains.

Recommended Dilution


WB: 1: 500 - 1: 2000

IHC: 1: 100 - 1: 300


ELISA: 1: 40000

Not yet tested in other applications.

Images

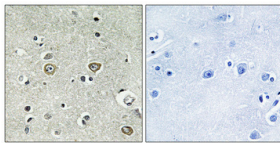

Immunohistochemical analysis of paraffin-embedded Human brain. Antibody was diluted at 1:100(4° overnight). High-pressure and temperature Tris-EDTA,pH8.0 was used for antigen retrieval. Negative contrl (right) obtained from antibody was pre-absorbed by immunogen peptide.

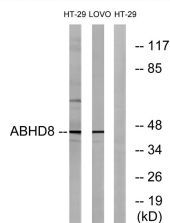

Western blot analysis of lysates from HT-29 and LOVO cells, using ABHD8 Antibody. The lane on the right is blocked with the synthesized peptide.

Storage

-20°C for one year