

FGF-BP2 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP03221
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from human FGFP2. AA range:170-219
Mol wt	24581
Species reactivity	Human
Clonality	Polyclonal
Recommended application	WB, ELISA
Concentration	1 mg/ml
Full name	FGF-BP2 Antibody
Synonyms	FGFBP2; KSP37; Fibroblast growth factor-binding protein 2; FGF-BP2; FGF-binding protein 2; FGFBP-2; 37 kDa killer-specific secretory protein; Ksp37; HBp17-related protein; HBp17-RP

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

FGFBP2 (fibroblast growth factor binding protein 2) encodes a member of the fibroblast growth factor binding protein family. The encoded protein is a serum protein that is selectively secreted by cytotoxic lymphocytes and may be involved in cytotoxic lymphocyte-mediated immunity. An increase in the amount of gene product may be associated with atopic asthma and mild extrinsic asthma.


Recommended Dilution

WB: 1: 500 - 1: 2000

ELISA: 1: 40000

Not yet tested in other applications.

Images

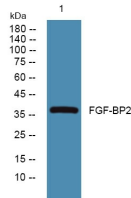

Western blot analysis of lysates from DU145 cells, primary antibody was diluted at 1:1000, 4° over night

Storage

-20°C for one year