

GABAA R Alpha1 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP03425
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from the Internal region of human GABRA1. AA range:61-110
Mol wt	50160
Species reactivity	Human, Mouse, Rat
Clonality	Polyclonal
Recommended application	WB, IHC-p, ELISA
Concentration	1 mg/ml
Full name	GABAA Ralpha1 Antibody
Synonyms	GABRA1; Gamma-aminobutyric acid receptor subunit alpha-1; GABA(A) receptor subunit alpha-1

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

GABRA1 encodes a gamma-aminobutyric acid (GABA) receptor. GABA is the major inhibitory neurotransmitter in the mammalian brain where it acts at GABA-A receptors, which are ligand-gated chloride channels. Chloride conductance of these channels can be modulated by agents such as benzodiazepines that bind to the GABA-A receptor. GABA-A receptors are pentameric, consisting of proteins from several subunit classes: alpha, beta, gamma, delta and rho. Mutations in this gene cause juvenile myoclonic epilepsy and childhood absence epilepsy type 4. Multiple transcript variants encoding the same protein have been identified for this gene.

Recommended Dilution


WB: 1: 500 - 1: 2000

IHC-p: 1: 100 - 1: 300


ELISA: 1: 10000

Not yet tested in other applications.


Images

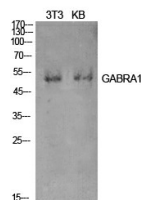


Western Blot analysis of NIH-3T3, KB cells using GABAA R α 1 Polyclonal Antibody. Secondary antibody was diluted at 1:20000

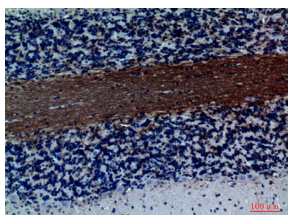

Immunohistochemical analysis of paraffin-embedded rat-brain, antibody was diluted at 1:100

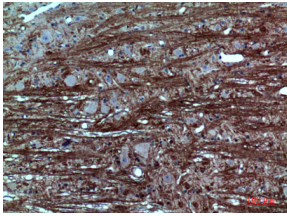

Immunohistochemical analysis of paraffin-embedded rat-brain, antibody was diluted at 1:100

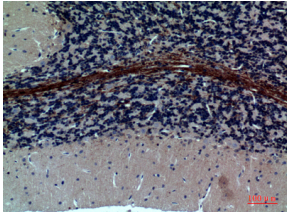

Immunohistochemical analysis of paraffin-embedded rat-brain, antibody was diluted at 1:100

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com