

## NMNA3 Rabbit Polyclonal Antibody

### Description

Product type	Primary Antibody
Code	BT-AP03833
Host	Rabbit
Isotype	IgG
Size	100ul, 50ul, 20ul
Immunogen	Synthesized peptide derived from human NMNA3
Mol wt	27720
Species reactivity	Human, Mouse
Clonality	Polyclonal
Recommended application	WB
Concentration	1 mg/ml
Full name	NMNA3
Synonyms	NMNA3

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

### Background


This gene encodes a member of the nicotinamide/nicotinic acid mononucleotide adenylyltransferase family. These enzymes use ATP to catalyze the synthesis of nicotinamide adenine dinucleotide or nicotinic acid adenine dinucleotide from nicotinamide mononucleotide or nicotinic acid mononucleotide respectively. The encoded protein is localized to mitochondria and may also play a neuroprotective role as a molecular chaperone. Alternatively spliced transcript variants encoding multiple isoforms have been observed for this gene.

### Recommended Dilution

WB: 1: 500 - 1: 2000

Not yet tested in other applications.

### Images

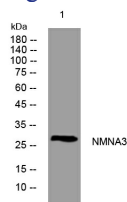

Western blot analysis of lysates from MCF-7 cells, primary antibody was diluted at 1:1000, 4°C overnight

### Storage

-20°C for 1 year