

HDAC1 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP03919
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from human HDAC1. AA range:387-436
Mol wt	55103
Species reactivity	Human, Mouse, Rat
Clonality	Polyclonal
Recommended application	WB, IHC-p, ELISA
Concentration	1 mg/ml
Full name	HDAC1 Antibody
Synonyms	HDAC1; RPD3L1; Histone deacetylase 1; HD1

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

Histone acetylation and deacetylation, catalyzed by multisubunit complexes, play a key role in the regulation of eukaryotic gene expression. Histone deacetylase 1 encoded by HDAC1 belongs to the histone deacetylase/acuc/apha family and is a component of the histone deacetylase complex. It also interacts with retinoblastoma tumor-suppressor protein and this complex is a key element in the control of cell proliferation and differentiation. Together with metastasis-associated protein-2, it deacetylates p53 and modulates its effect on cell growth and apoptosis.

Recommended Dilution


WB: 1: 500 - 1: 2000

IHC: 1: 100 - 1: 300


ELISA: 1: 5000

Not yet tested in other applications.


Images

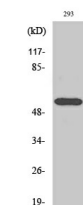


Western Blot analysis of various cells using HDAC1 Polyclonal Antibody diluted at 1:500

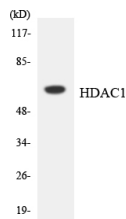

Western blot analysis of the lysates from COLO205 cells using HDAC1 antibody.

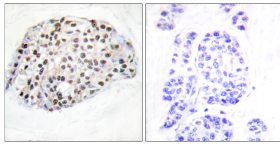

Immunohistochemistry analysis of paraffin-embedded human breast carcinoma tissue, using HDAC1 Antibody. The picture on the right is blocked with the synthesized peptide.

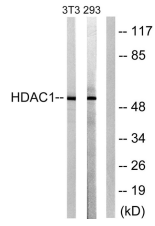

Western blot analysis of lysates from NIH/3T3 cells, using HDAC1 Antibody. The lane on the right is blocked with the synthesized peptide.

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com