

S2545 Rabbit Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP04593
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	Synthesized peptide derived from human S2545
Mol wt	31680
Species reactivity	Human, Mouse
Clonality	Polyclonal
Recommended application	WB
Concentration	1 mg/ml
Full name	S2545
Synonyms	S2545

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background


SLC25A45 belongs to the SLC25 family of mitochondrial carrier proteins (Haitina et al. 2006)

Recommended Dilution

WB: 1: 500 - 1: 2000

Not yet tested in other applications.

Images

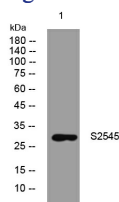

Western blot analysis of lysates from 3T3 cells, primary antibody was diluted at 1:1000, 4°C overnight

Storage

-20°C for 1 year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com