

MRP-L13 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP05571
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from human MRPL13. AA range:121-170
Mol wt	20692
Species reactivity	Human, Mouse, Rat
Clonality	Polyclonal
Recommended application	WB, IHC-p, ELISA
Concentration	1 mg/ml
Full name	MRP-L13 Antibody
Synonyms	MRPL13; 39S ribosomal protein L13; mitochondrial; L13mt; MRP-L13

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

Mammalian mitochondrial ribosomal protein L 13 (39S ribosomal protein L13 mitochondrial) is encoded by nuclear gene MRPL13 and help in protein synthesis within the mitochondrion. Mitochondrial ribosomes (mitoribosomes) consist of a small 28S subunit and a large 39S subunit. They have an estimated 75% protein to rRNA composition compared to prokaryotic ribosomes, where this ratio is reversed. Another difference between mammalian mitoribosomes and prokaryotic ribosomes is that the latter contain a 5S rRNA. Among different species, the proteins comprising the mitoribosome differ greatly in sequence, and sometimes in biochemical properties, which prevents easy recognition by sequence homology.

Recommended Dilution

WB: 1: 500 - 1: 2000

IHC: 1: 100 - 1: 300


ELISA: 1: 40000

Not yet tested in other applications.


Images



Western Blot analysis of various cells using MRP-L13 Polyclonal Antibody diluted at 1:1000

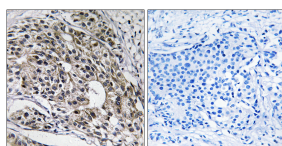

Immunohistochemistry analysis of paraffin-embedded human breast carcinoma tissue, using MRPL13 Antibody. The picture on the right is blocked with the synthesized peptide.

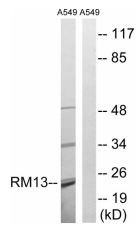

Western blot analysis of lysates from A549 cells, using MRPL13 Antibody. The lane on the right is blocked with the synthesized peptide.

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com