

Olfactory receptor 5AP2 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP06499
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from human OR5AP2. AA range:267-316
Mol wt	35508
Species reactivity	Human
Clonality	Polyclonal
Recommended application	WB, IF, ELISA
Concentration	1 mg/ml
Full name	Olfactory receptor 5AP2 Antibody
Synonyms	OR5AP2; Olfactory receptor 5AP2

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

Olfactory receptors interact with odorant molecules in the nose, to initiate a neuronal response that triggers the perception of a smell. The olfactory receptor proteins are members of a large family of G-protein-coupled receptors (GPCR) arising from single coding-exon genes. Olfactory receptors share a 7-transmembrane domain structure with many neurotransmitter and hormone receptors and are responsible for the recognition and G protein-mediated transduction of odorant signals. The olfactory receptor gene family is the largest in the genome. The nomenclature assigned to the olfactory receptor genes and proteins for this organism is independent of other organisms.

Recommended Dilution


WB: 1: 500 - 1: 2000

IF: 1: 200 - 1: 1000


ELISA: 1: 5000

Not yet tested in other applications.

Images

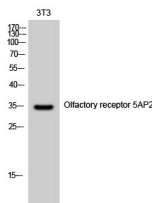

Western Blot analysis of 3T3 cells using Olfactory receptor 5AP2 Polyclonal Antibody

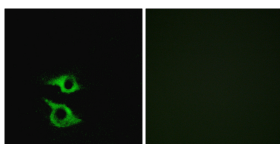

Immunofluorescence analysis of HeLa cells, using OR5AP2 Antibody. The picture on the right is blocked with the synthesized peptide.

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com