

Olfactory receptor 8J1 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP06553
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from human OR8J1. AA range:233-282
Mol wt	35355
Species reactivity	Human
Clonality	Polyclonal
Recommended application	WB, IF, ELISA
Concentration	1 mg/ml
Full name	Olfactory receptor 8J1 Antibody
Synonyms	OR8J1; Olfactory receptor 8J1; Olfactory receptor OR11-183

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

Olfactory receptors interact with odorant molecules in the nose, to initiate a neuronal response that triggers the perception of a smell. The olfactory receptor proteins are members of a large family of G-protein-coupled receptors (GPCR) arising from single coding-exon genes. Olfactory receptors share a 7-transmembrane domain structure with many neurotransmitter and hormone receptors and are responsible for the recognition and G protein-mediated transduction of odorant signals. The olfactory receptor gene family is the largest in the genome. The nomenclature assigned to the olfactory receptor genes and proteins for this organism is independent of other organisms.

Recommended Dilution


WB: 1: 500 - 1: 2000

IF: 1: 200 - 1: 1000


ELISA: 1: 10000

Not yet tested in other applications.


Images

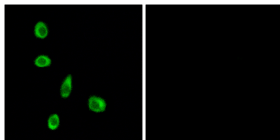

Immunofluorescence analysis of MCF7 cells, using OR8J1 Antibody. The picture on the right is blocked with the synthesized peptide.

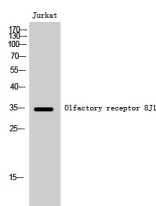

Western Blot analysis of Jurkat cells using Olfactory receptor 8J1 Polyclonal Antibody

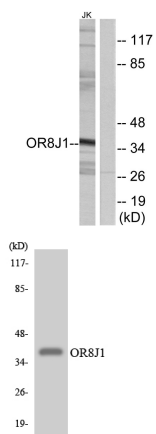

Western blot analysis of lysates from Jurkat cells, using OR8J1 Antibody. The lane on the right is blocked with the synthesized peptide.

Western blot analysis of the lysates from HepG2 cells using OR8J1 antibody.

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com