

SHIP-2 Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP08296
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	The antiserum was produced against synthesized peptide derived from the Internal region of human INPPL1. AA range:351-400
Mol wt	138585
Species reactivity	Human, Mouse, Rat
Clonality	Polyclonal
Recommended application	WB, IHC-p, ELISA
Concentration	1 mg/ml
Full name	SHIP-2 Antibody
Synonyms	INPPL1; SHIP2; Phosphatidylinositol 3,4,5-trisphosphate 5-phosphatase 2; Inositol polyphosphate phosphatase-like protein 1; INPPL-1; Protein 51C; SH2 domain-containing inositol 5'-phosphatase 2; SH2 d

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

Phosphatidylinositol 3,4,5-trisphosphate 5-phosphatase 2 encoded by INPPL1 is an SH2-containing 5'-inositol phosphatase that is involved in the regulation of insulin function. The encoded protein also plays a role in the regulation of epidermal growth factor receptor turnover and actin remodelling. Additionally, this gene supports metastatic growth in breast cancer and is a valuable biomarker for breast cancer.

Recommended Dilution


WB: 1: 500 - 1: 2000

IHC-p: 1: 100 - 300


ELISA: 1: 20000

Not yet tested in other applications.


Images

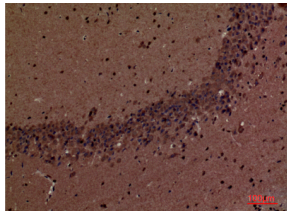

Immunohistochemical analysis of paraffin-embedded mouse-brain, antibody was diluted at 1:100

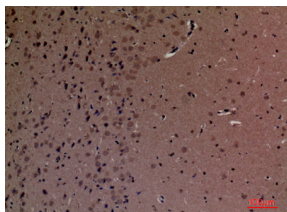

Immunohistochemical analysis of paraffin-embedded mouse-brain, antibody was diluted at 1:100

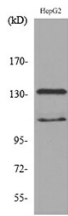

Western blot analysis of lysate from HepG2 cells, using INPPL1 Antibody.

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com