

## SLAM Polyclonal Antibody

### Description

|                                |  |
|--------------------------------|--|
| <b>Product type</b>            | Primary Antibody   |
| <b>Code</b>                    | BT-AP08348   |
| <b>Host</b>                    | Rabbit   |
| <b>Isotype</b>                 | IgG  |
| <b>Size</b>                    | 20ul, 50ul, 100ul  |
| <b>Immunogen</b>               | The antiserum was produced against synthesized peptide derived from the Internal region of human SLAMF1. AA range:81-130 |
| <b>Mol wt</b>                  | 37231  |
| <b>Species reactivity</b>      | Human  |
| <b>Clonality</b>               | Polyclonal   |
| <b>Recommended application</b> | WB, IHC-p, ELISA   |
| <b>Concentration</b>           | 1 mg/ml  |
| <b>Full name</b>               | SLAM Antibody  |
| <b>Synonyms</b>                | SLAMF1; SLAM; Signaling lymphocytic activation molecule; CDw150; IPO-3; CD150  |

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

### Background

High-affinity self-ligand important in bidirectional T-cell to B-cell stimulation. SLAM-induced signal-transduction events in T-lymphocytes are different from those in B-cells. Two modes of SLAM signaling are likely to exist: one in which the inhibitor SH2D1A acts as a negative regulator and another in which protein-tyrosine phosphatase 2C (PTPN11)-dependent signal transduction operates.

### Recommended Dilution


WB: 1: 500 - 1: 2000

IHC-p: 1: 100 - 300


ELISA: 1: 20000

Not yet tested in other applications.


### Images

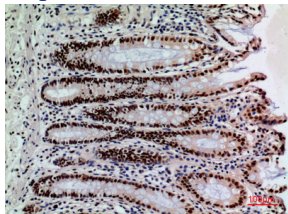

Immunohistochemical analysis of paraffin-embedded human-colon, antibody was diluted at 1:100

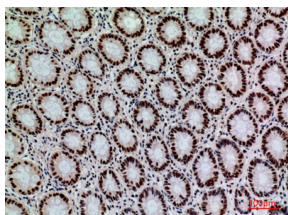


Immunohistochemical analysis of paraffin-embedded human-colon, antibody was diluted at 1:100

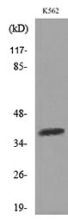

Western blot analysis of lysate from K562 cells, using SLAMF1 Antibody.



Western Blot analysis of K562 cells using SLAM Polyclonal Antibody. Secondary antibody was diluted at 1:20000

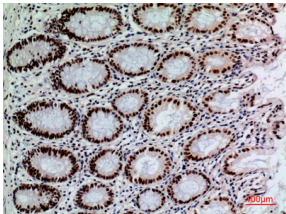

Immunohistochemical analysis of paraffin-embedded human-colon, antibody was diluted at 1:100

### Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: [save@bt-laboratory.com](mailto:save@bt-laboratory.com) | [www.bt-laboratory.com](http://www.bt-laboratory.com)