

DAD1 Rabbit Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP08903
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	Synthesized peptide derived from human DAD1
Mol wt	12430
Species reactivity	Human, Mouse, Rat
Clonality	Polyclonal
Recommended application	WB
Concentration	1 mg/ml
Full name	DAD1
Synonyms	DAD1

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background


DAD1| the defender against apoptotic cell death| was initially identified as a negative regulator of programmed cell death in the temperature sensitive tsBN7 cell line. The DAD1 protein disappeared in temperature-sensitive cells following a shift to the nonpermissive temperature| suggesting that loss of the DAD1 protein triggered apoptosis. DAD1 is believed to be a tightly associated subunit of oligosaccharyltransferase both in the intact membrane and in the purified enzyme| thus reflecting the essential nature of N-linked glycosylation in eukaryotes.

Recommended Dilution

WB: 1: 500 - 1: 2000

Not yet tested in other applications.

Images

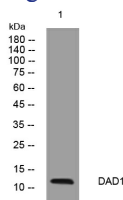

Western blot analysis of lysates from MDA-MB cells, primary antibody was diluted at 1:1000, 4°C overnight

Storage

-20°C for 1 year