

NDUAB Polyclonal Antibody

Description

Product type	Primary Antibody
Code	BT-AP11756
Host	Rabbit
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	Synthesized peptide derived from human protein . at AA range: 30-110
Mol wt	N/A
Species reactivity	Human, Rat, Mouse
Clonality	Polyclonal
Recommended application	WB, ELISA
Concentration	1 mg/ml
Full name	NADH dehydrogenase [ubiquinone] 1 alpha subcomplex subunit 11
Synonyms	NADH dehydrogenase [ubiquinone] 1 alpha subcomplex subunit 11 ;Complex I-B14.7;CI-B14.7;NADH-ubiquinone oxidoreductase subunit B14.7

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

This gene encodes a subunit of the membrane-bound mitochondrial complex I. Complex I is composed of numerous subunits and functions as the NADH-ubiquinol reductase of the mitochondrial electron transport chain. Mutations in this gene are associated with severe mitochondrial complex I deficiency. Alternate splicing results in multiple transcript variants.


Recommended Dilution

WB: 1: 500 - 1: 2000

ELISA: 1: 5000 - 1: 20000

Not yet tested in other applications.

Images

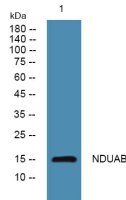

Western blot analysis of lysates from Jarkat cells, primary antibody was diluted at 1:1000, 4°C overnight

Storage

-20°C for 1 year