

Survivin Monoclonal Antibody(Mix)

Description

Product type	Primary Antibody
Code	BT-MCA0074
Host	Mouse
Isotype	IgG
Size	20ul, 50ul, 100ul
Immunogen	Recombinant Protein of Survivin
Mol wt	N/A
Species reactivity	Human,Rat
Clonality	Monoclonal
Recommended application	WB, IHC-p, IF
Concentration	1 mg/ml
Full name	Baculoviral IAP repeat-containing protein 5
Synonyms	BIRC5; API4; IAP4; Baculoviral IAP repeat-containing protein 5; Apoptosis inhibitor 4; Apoptosis inhibitor survivin

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

This gene is a member of the inhibitor of apoptosis (IAP) gene family, which encode negative regulatory proteins that prevent apoptotic cell death. IAP family members usually contain multiple baculovirus IAP repeat (BIR) domains, but this gene encodes proteins with only a single BIR domain. The encoded proteins also lack a C-terminus RING finger domain. Gene expression is high during fetal development and in most tumors, yet low in adult tissues. Alternatively spliced transcript variants encoding distinct isoforms have been found for this gene.


Recommended Dilution

IHC: 1:200-500


WB: 1:1000-2000

Not yet tested in other applications.


Images

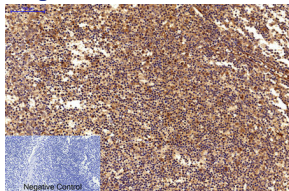


Immunohistochemical analysis of paraffin-embedded Human-Tonsil tissue. 1.Survivin Monoclonal antibody(Mix) was diluted at 1:200(4°C,overnight). 2.Sodium citrate pH 6.0 was used for antibody retrieval(>98°C,20min). 3.Secondary antibody was diluted at 1:200(room temperature, 30min). Negative control was used by secondary antibody only.

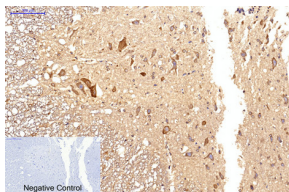


Immunohistochemical analysis of paraffin-embedded Rat-spinal-cord tissue. 1.Survivin Monoclonal antibody(Mix) was diluted at 1:200(4°C,overnight). 2.Sodium citrate pH 6.0 was used for antibody retrieval(>98°C,20min). 3.Secondary antibody was diluted at 1:200(room temperature, 30min). Negative control was used by secondary antibody only.

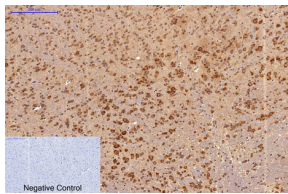


Immunohistochemical analysis of paraffin-embedded Mouse-brain tissue. 1.Survivin Monoclonal antibody(Mix) was diluted at 1:200(4°C,overnight). 2.Sodium citrate pH 6.0 was used for antibody retrieval(>98°C,20min). 3.Secondary antibody was diluted at 1:200(room temperature, 30min). Negative control was used by secondary antibody only.

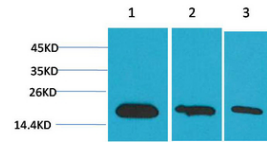


Western blot analysis of 1) HeLa, 2) 293, 3) PC12 using Survivin Monoclonal antibody.

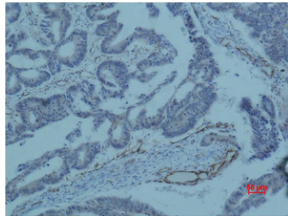


Immunohistochemical analysis of paraffin-embedded human colon carcinoma using Survivin Monoclonal antibody.

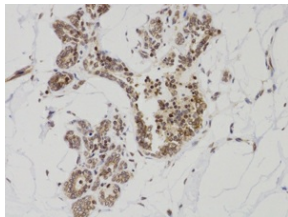


Immunohistochemical analysis of paraffin-embedded Human breast cancer. 1. Using Survivin Mouse Monoclonal antibody (YM3419) diluted at 1:200 (4°C overnight).2.High-pressure and temperature Citric acid, pH6.0 was used for antigen retrieval. 3.Secondary antibody was diluted at 1:200(room temperature, 50min). Picture was kindly provided by our customer from Tianjin Medical University Cancer Institute and Hospital

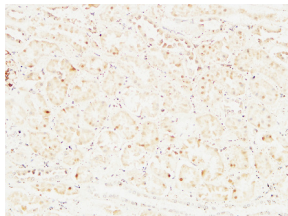

Immunohistochemical analysis of paraffin-embedded Human kidney.1.Antibody was diluted at 1:200(4°C overnight). 2.High-pressure and temperature EDTA, pH8.0 was used for antigen retrieval. 3.Secondary antibody was diluted at 1:200(room temperature, 30min).

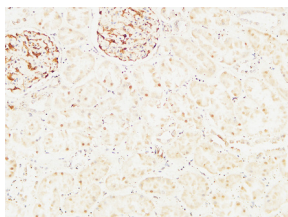

Immunohistochemical analysis of paraffin-embedded Human kidney.1.Antibody was diluted at 1:200(4°C overnight). 2.High-pressure and temperature EDTA, pH8.0 was used for antigen retrieval. 3.Secondary antibody was diluted at 1:200(room temperature, 30min).

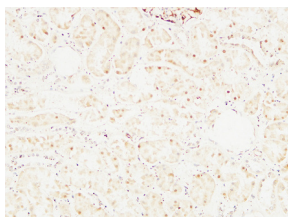

Immunohistochemical analysis of paraffin-embedded Human kidney.1.Antibody was diluted at 1:200(4°C overnight). 2.High-pressure and temperature EDTA, pH8.0 was used for antigen retrieval. 3.Secondary antibody was diluted at 1:200(room temperature, 30min).

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com