

NGFR p75 Monoclonal Antibody

Description

| | |
|--------------------------------|---|
| Product type | Primary Antibody |
| Code | BT-MCA0114 |
| Host | Mouse |
| Isotype | IgG |
| Size | 50ul, 100ul |
| Immunogen | Purified recombinant fragment of human NGFR p75 expressed in E. Coli. |
| Mol wt | N/A |
| Species reactivity | Human |
| Clonality | Monoclonal |
| Recommended application | WB, IF, ICC, FCM, ELISA |
| Concentration | 1 mg/ml |
| Full name | Tumor necrosis factor receptor superfamily member 16 |
| Synonyms | NGFR; TNFRSF16; Tumor necrosis factor receptor superfamily member 16; Gp80-LNGFR; Low affinity neurotrophin receptor p75NTR; Low-affinity nerve growth factor receptor; NGF receptor; p75 ICD; CD antigen CD271 |

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

Nerve growth factor receptor contains an extracellular domain containing four 40-amino acid repeats with 6 cysteine residues at conserved positions followed by a serine/threonine-rich region, a single transmembrane domain, and a 155-amino acid cytoplasmic domain. The cysteine-rich region contains the nerve growth factor binding domain.

Recommended Dilution

ELISA: 1:10000


FC: 1:200 - 1:400

IF: 1:200 - 1:1000


WB: 1:500 - 1:2000

Not yet tested in other applications.


Images

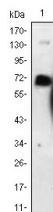

Western Blot analysis using NGFR p75 Monoclonal antibody against NGFR-hlgGfC transfected HEK293 cell lysate.

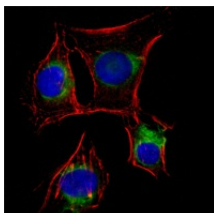

Immunofluorescence analysis of EC cells using NGFR p75 Monoclonal antibody (green). Red: Actin filaments have been labeled with DY-554 phalloidin. Blue: DRAQ5 fluorescent DNA dye.

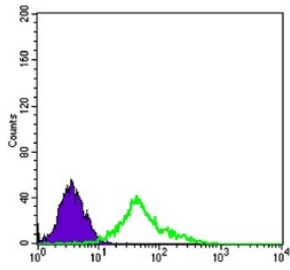

Flow cytometric analysis of EC cells using NGFR p75 Monoclonal antibody (green) and negative control (purple).

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com