

CD1A Monoclonal Antibody

Description

Product type	Primary Antibody
Code	BT-MCA0283
Host	Mouse
Isotype	IgG
Size	50ul, 100ul
Immunogen	Purified recombinant fragment of human CD1A expressed in E. Coli.
Mol wt	N/A
Species reactivity	Human
Clonality	Monoclonal
Recommended application	WB, FCM, ELISA
Concentration	1 mg/ml
Full name	T-cell surface glycoprotein CD1a
Synonyms	CD1A; T-cell surface glycoprotein CD1a; T-cell surface antigen T6; Leu-6; hTa1 thymocyte antigen; CD antigen CD1a

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

This gene encodes a member of the CD1 family of transmembrane glycoproteins, which are structurally related to the major histocompatibility complex (MHC) proteins and form heterodimers with beta-2-microglobulin. The CD1 proteins mediate the presentation of primarily lipid and glycolipid antigens of self or microbial origin to T cells. The human genome contains five CD1 family genes organized in a cluster on chromosome 1. The CD1 family members are thought to differ in their cellular localization and specificity for particular lipid ligands. The protein encoded by this gene localizes to the plasma membrane and to recycling vesicles of the early endocytic system. Alternative splicing results in multiple transcript variants.

Recommended Dilution


ELISA: 1:10000

FC: 1:200 - 1:400


WB: 1:500 - 1:2000

Not yet tested in other applications.

Images

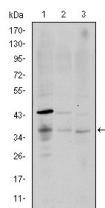

Western Blot analysis using CD1A Monoclonal antibody against K562 (1) RAJI (2) and MOLT4 (3) cell lysate.

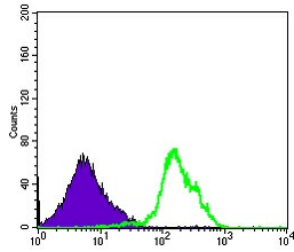

Flow cytometric analysis of RAJI cells using CD1A Monoclonal antibody (green) and negative control (purple).

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com