

Cyclin D3 Monoclonal Antibody

Description

Product type	Primary Antibody
Code	BT-MCA0397
Host	Mouse
Isotype	IgG
Size	50ul, 100ul
Immunogen	Purified recombinant fragment of human Cyclin D3 expressed in E. Coli.
Mol wt	N/A
Species reactivity	Human,Mouse
Clonality	Monoclonal
Recommended application	WB, ELISA
Concentration	1 mg/ml
Full name	G1,S-specific cyclin-D3
Synonyms	CCND3; G1; S-specific cyclin-D3

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

The protein encoded by this gene belongs to the highly conserved cyclin family, whose members are characterized by a dramatic periodicity in protein abundance through the cell cycle. Cyclins function as regulators of CDK kinases. Different cyclins exhibit distinct expression and degradation patterns which contribute to the temporal coordination of each mitotic event. This cyclin forms a complex with and functions as a regulatory subunit of CDK4 or CDK6, whose activity is required for cell cycle G1/S transition. This protein has been shown to interact with and be involved in the phosphorylation of tumor suppressor protein Rb. The CDK4 activity associated with this cyclin was reported to be necessary for cell cycle progression through G2 phase into mitosis after UV radiation. Several transcript variants encoding different isoforms have been found for this gene.


Recommended Dilution

ELISA: 1:10000

WB: 1:500 - 1:2000

Not yet tested in other applications.

Images

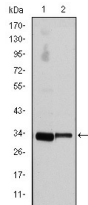

Western Blot analysis using Cyclin D3 Monoclonal antibody against NIH/3T3 (1) and Jurkat (2) cell lysate.

Storage

-20°C for one year