

HDAC3 Monoclonal Antibody

Description

Product type	Primary Antibody
Code	BT-MCA0667
Host	Mouse
Isotype	IgG
Size	50ul, 100ul
Immunogen	Purified recombinant fragment of human HDAC3 (aa224-428) expressed in E. Coli.
Mol wt	N/A
Species reactivity	Human
Clonality	Monoclonal
Recommended application	WB, IF, ICC, ELISA
Concentration	1 mg/ml
Full name	Histone deacetylase 3
Synonyms	HDAC3; Histone deacetylase 3; HD3; RPD3-2; SMAP45

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

Histones play a critical role in transcriptional regulation, cell cycle progression, and developmental events. Histone acetylation/deacetylation alters chromosome structure and affects transcription factor access to DNA. The protein encoded by this gene belongs to the histone deacetylase/acuc/apha family. It has histone deacetylase activity and represses transcription when tethered to a promoter. It may participate in the regulation of transcription through its binding with the zinc-finger transcription factor YY1. This protein can also down-regulate p53 function and thus modulate cell growth and apoptosis. This gene is regarded as a potential tumor suppressor gene.

Recommended Dilution


ELISA: 1:10000

IF: 1:200 - 1:1000


WB: 1:500 - 1:2000

Not yet tested in other applications.

Images

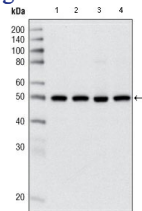

Western Blot analysis using HDAC3 Monoclonal antibody against HeLa (1) NIH/3T3 (2) C6 (3) and COS (4) cell lysate.

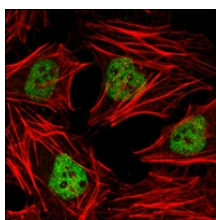

Confocal immunofluorescence analysis of HeLa cells using HDAC3 Monoclonal antibody (green). Red: Actin filaments have been labeled with DY-554 phalloidin.

Storage

-20°C for one year

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com