



Competitive ELISA Kit

User Instruction

Sheep Progesterone, PROG ELISA Kit

Cat.No: EA0004Sh

CONTENTS

Precision.....	4
Intended Use.....	4
Assay Principle.....	4
Reagents Provided.....	5
Materials Required but Not Supplied.....	5
Precautions.....	6
Specimen Collection.....	6
Reagents Preparation.....	7
Assay Procedure.....	8
Summary.....	9
Calculation of Results.....	9
Typical Data.....	10
Troubleshooting.....	10
Data Template.....	12

Precision

Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays) Three samples of known concentration were tested in separate assays to assess inter-assay precision.

CV (%) = SD/mean x 100

Intra-assay: CV < 10%

Inter-assay: CV< 12%

Standard curve range: 0.2 - 60ng/ml

Sensitivity: 0.097ng/ml

Size: 96 Test

Storage: Store at -20°C for one year. Or store at 2-8°C for 6 months. If individual reagents are opened it is recommended that the kit be used within 1 month. Avoid repeated thaw cycles.


Intended Use

This Competitive ELISA kit is for the quantification of Sheep Progesterone (also known as PROG) in serum, plasma, cell culture supernatants, urine, tissue homogenates or other biological fluids.

Note: This product is for research use only, not for use in diagnosis procedures. It's highly recommended to read this instruction entirely before use.

Assay Principle

This kit is an Enzyme-Linked Immunosorbent Assay (ELISA). Add samples to the pre-coated plate. Then add biotinylated antigen. The antigens in the samples compete with the biotinylated antigen to bind to the capture antibody and incubate. Unbound antigen is washed away during a washing step. An avidin-HRP is then added and then incubate. Unbound avidin-HRP is washed away during a washing step. TMB Substrate is then added and color develops. The reaction is stopped by addition of acidic stop solution and color changes into yellow that can be measured at 450 nm. The intensity of the color developed is inversely proportional to the concentration of PROG in the sample. The concentration of PROG in the sample is then determined by comparing the O.D. of the samples to the standard curve.

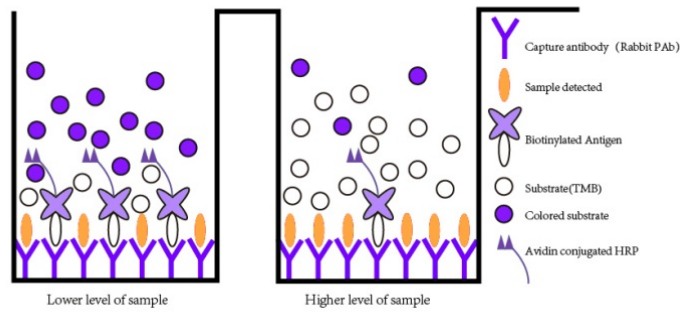

Reagent Provided

Components	Quantity
Pre-coated plate	12 * 8 well strips x 1
Sheep PROG Standard, lyophilized	2 vials
Standard/Sample diluent	6ml x1 vial
Biotinylated antigen, lyophilized	1 vial
Avidin-HRP concentrate	100ul x 1 vial
Biotinylated antigen diluent	6ml x1
Avidin HRP diluent	5.9ml x 1 vials
Substrate solution A	6ml x1
Substrate solution B	6ml x1
Stop solution	6ml x1
Wash buffer (25x)	20ml x 1 vial
Plate sealer	2 pics
User instruction	1

Material Required but Not Supplied

- 37°C±0.5°C incubator

- Precision pipette and disposable tip
- Deionized or distilled water
- Clean tubes
- Absorbent paper
- Microplate reader with $450 \pm 10\text{nm}$ wavelength filter

Precautions

- Prior to running the assay, the kit and sample should be warmed naturally to room temperature 30 minutes.
- Once the desired number of strips has been removed, immediately reseal the bag to protect the remain from deterioration. Cover all reagents when not in use.
- Make sure pipetting order and rate of addition from well-to-well when pipetting reagents.
- There are two vials of standard in the kit for users, please cover the other unused vial and keep refrigerated.
- Do not allow wells to become dry during the assay procedure.
- This instruction should be strictly followed in the experiment.
- Pipette tips and plate sealer in hand should be clean and disposable to avoid cross-contamination.
- Avoid using the reagents from different batches together.
- Substrate solution B is sensitive to light, don't expose substrate solution B to light for a long time.
- Stop solution contains acid. Please wear eye, hand and skin protection when using this material. Avoid contact of skin or mucous membranes with kit reagent.
- The kit should not be used beyond the expiration date.

Specimen Collection

Serum Allow serum to clot for 10-20 minutes at room temperature. Centrifuge at 2000-3000 RPM for 20 minutes. Collect the supernatant without sediment.

Plasma Collect plasma using EDTA or heparin as an anticoagulant. After mix 10-20 minutes, centrifuge samples for 20 minutes at 2000-3000 RPM. Collect the supernatant without sediment.

Cell culture supernatant Collect by sterile tubes. When detecting secrete components, centrifuge at 2000-3000 RPM for 20 minutes. Collect the supernatants. When detecting the components in the cell, use PBS (pH 7.2-7.4) to dilute cell suspension, the cell concentration of approximately 1 million/ml. Damage cells through repeated freeze-thaw cycles to let out the inside components. Centrifuge at 2000-3000 RPM for 20 minutes. Collect the supernatant without sediment.

Tissue Rinse tissues in ice-cold PBS (pH 7.4) to remove excess blood thoroughly and weigh before homogenization. Mince tissues and homogenize them in PBS (tissue weight (g): PBS (mL) volume=1:9) with a glass homogenizer on ice. To further break down the cells, you can sonicate the suspension with an ultrasonic cell disrupter or subject it to freeze-thaw cycles. The homogenates are then centrifuged for 15 minutes at 12,000 RPM at 4°C to get the supernatant. Avoid freeze/thaw cycles.


Urine/Ascites/Cerebrospinal fluid Collect by sterile tube. Centrifuge at 2000-3000 RPM for 20 minutes. Collect the supernatant without sediment.

Note

- The supplier is only responsible for the kit itself, but not for the samples consumed during the assay. The user should calculate the possible amount of the sample used in the whole test
- Samples to be used within 5 days should be stored at 2-8°C. Samples should be aliquoted or must be stored at -20°C within 1 month or -80°C within 3 months. Avoid repeated freeze thaw cycles.
- Samples should be brought to room temperature before starting the assay.
- Samples containing NaN₃ can't be tested as it inhibits the activity of Horse Radish Peroxidase (HRP).
- Collect the supernatants carefully. When sediments occurred during storage, centrifugation should be performed again.
- Hemolysis can greatly impact the validity of test results. Take care to minimize hemolysis.

Reagent Preparation

- All reagents should be brought to room temperature before use.
- **Standard** Reconstitute one vial of standard with 150ul of Standard/Sample diluent to generate a 64ng/ml standard stock solution which should be used within 24 hours. Allow the standard to sit for 15 minutes with gentle agitation prior to making dilutions. Prepare duplicate or triplicate standard points by serially diluting the standard stock solution 1: 2 with diluent to produce 32ng/ml, 16ng/ml, 8ng/ml, 4ng/ml and 2ng/ml solutions. Add Standard/Sample diluent only as the zero standard (0ng/ml).

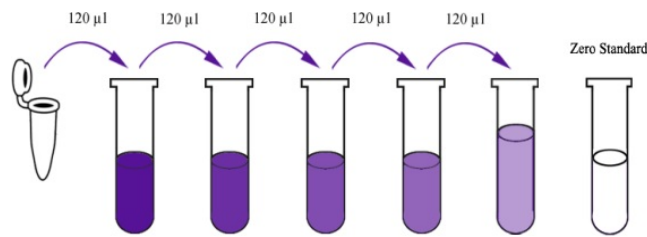

Stock solution	64ng/ml	32ng/ml	16ng/ml	8ng/ml	4ng/ml	2ng/ml
----------------	---------	---------	---------	--------	--------	--------

- **Biotinylated Antigen** Briefly centrifuge the biotinylated antigen vial then add 1ml Biotinylated antigen diluent to mix well. And then pipette all this solution back into the Biotinylated antigen diluent vial to mix well and generate a 6ml stock solution. Allow to sit for 10 minutes with gentle agitation prior to making dilutions.
- **Avidin-HRP Concentrate** Briefly low- speed centrifuge the avidin-HRP Concentrates solution and then pipette all avidin-HRP into the Avidin HRP Diluent vial. Mix well to generate a 6ml stock solution. Allow to sit for 10 minutes with gentle agitation prior to making dilutions.
- **Wash Buffer Concentrate 25x** Dilute 20ml of concentrated wash buffer with 480ml double distilled water to prepare 500 ml of wash buffer. If crystals have formed in the concentrate, warm it in a 40°C water bath and mix it gently until the crystals have completely dissolved.

Assay Procedure

1. Prepare all reagents, standard solutions and samples as instructed. Bring all reagents to room temperature before use. The assay should be performed at room temperature.
2. Determine the number of strips required for the assay. Insert the strips in the frames for use. The unused strips should be stored at 4°C for up to one month.
3. Set a blank well without any solution.
4. Add 50ul negative control to each of the negative control wells and 50ul positive control to each of the positive control wells. Add 40ul sample diluent and then add 10ul sample to the sample well, mix well.
5. Remove the sealer and wash the plate 5 times with wash buffer. Soak wells with 300ul wash buffer for 30 seconds to 1 minute for each wash. For automated washing, aspirate all wells and wash 5 times with wash buffer. Blot the plate onto paper towels or other absorbent material.
6. Add 50ul HRP to each well (except blank well). Cover with a plate sealer, and incubate for 30 minutes at 37°C.
7. Remove the sealer and wash as described above.
8. Add 50ul substrate solution A to each well and then add 50ul substrate solution B to each well. Mix well. Incubate plate covered with a new sealer for 10 minutes at 37°C in the dark.
9. Add 50ul stop solution to each well, the blue color will change into yellow immediately.
10. Determine the optical density (OD value) of each well immediately using a microplate reader set to 450 nm within 15 minutes after adding the stop solution.

Summary

Prepare all reagents, samples and standards.

Add samples, standard and biotinylated antigen. Incubate for 60 minutes at 37°C.

Aspirate and wash 5 times.

Add avidin-HRP and incubate for 60 minutes at 37°C.

Aspirate and wash 5 times.

Add substrate solution A and substrate solution B and incubate in the dark for 10 minutes at 37°C.

Add stop solution.

Read the OD value within 10 minutes at 450nm.

Add 50ul stop solution.


Calculation of Result

Average the duplicate readings for each standard, control, and sample. Create a standard curve by plotting the mean absorbance for each standard on the Y-axis against the target antigen concentration on the X-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the target antigen concentration on the X axis versus the O.D. of the standards on the Y axis and the best fit line can be determined by regression analysis. The linear equation ($X = Y + \text{Calibration Value}$) can be used to calculate the standard curve where X is the log of the concentration of the standard and Y is the OD value of the standard. If samples have been diluted (2-5 times is recommended), the concentration read from the standard curve must be multiplied by the dilution factor.

Concentration	Blank	S5	S4	S3	S2	S1	S0
OD value	VB	V5	V4	V3	V2	V1	V0
Calibration value	VB-VB	V5-VB	V4-VB	V3-VB	V2-VB	V1-VB	V0-VB

Typical Data

This standard curve is only for demonstration purposes. A standard curve should be generated with each assay.

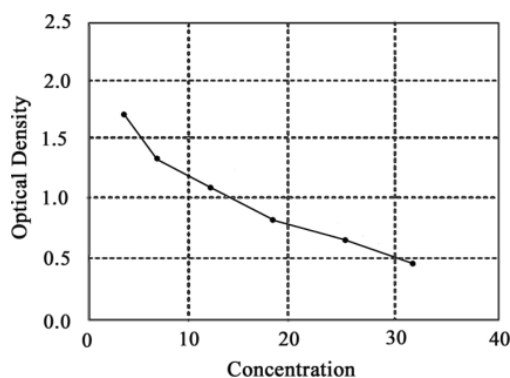

Troubleshooting

High Background possible case	Solution
Improper washing	Increasing duration of soaking steps
Incorrect incubation temperature	Incubate at 37°C
Incubation time too long	Reduce incubation time
Substrate exposed to light prior to use	Keep substrate in a dark place
Substrate was contaminated	Replace substrate. Substrate should be clean and avoid crossed contamination by using the sealer
Contaminated wash buffer	Use a clean buffers and sterile filter


Weak or No Signal possible case	Solution
A reagent or a step of the procedure omitted by mistake	Check protocol and follow steps carefully

Antibody are not enough	Increase the concentration of the antibody
Improper washing	Increasing duration of soaking steps
Reagent are contaminated	Use new one
Pipette are not clean	Pipette should be clean
HRP was not added	Add HRP according to the instruction
Sample contains sodium azide	Don't prepare samples with sodium azide
Wrong incubation time or temperature	Check and follow protocol. Place plates in an incubator during incubation periods (set to 37°C).

Poor Precision possible case	Solution
Pipetting error	Dispense quickly and identically into the side of each well. Use calibrated pipettes.
Incomplete washing	Make sure wells are washed adequately by filling the wells with proper amount of wash buffer.
Unclean wells	Inspect wells and remove debris prior to use. Wipe the outer bottom of plate clean to remove any debris or fingerprints prior to reading.

Poor standard curve possible case	Solution
Incorrect preparation of standard	Reconstitute standard as suggested on data sheet.
Capture Antibody did not bind	Use BT LAB ELISA plate in the kit
Inefficient washing	Be sure wash apparatus is working properly (i.e. distributing even volumes into each well). Be sure wells are empty after aspiration, yet be sure to fill wells in a timely manner.
Pipetting error	Dispense quickly and identically into the side of each well. Use calibrated pipettes.
Incorrect storage of components	Store all components as recommended on data sheet. Do not allow reconstituted reagents to stay at room temperature for excess time.

Assay Principle

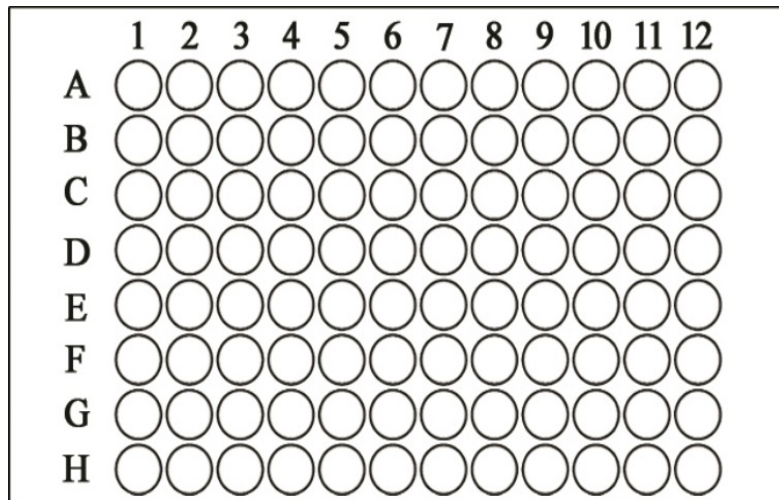

If you have any question on the order please contact us via: order@bt-laboratory.com;
 technical assistance please contact us via: support@bt-laboratory.com
 More product visit www.bt-laboratory.com